

# TR 900A

## Height adjustable Autofill bathtubs

Hygiene Equipment  
Designed for People

### Quick Guide for TR 900A, -CA, -CWA



#### Before and during bathing

Anyone operating the Bath/Lifter must fully understand all proceedings before - during and after the entire bathing cycle (see operating manuals).

#### Preparing Tub

- Press and hold control button (1) to lower the bathtub to its lowest position.
- Close the drain plug (9) by turning the overflow drainage knob.
- Set foot rest (7) into position appropriate for the patient legs. Footrest is not used with TR's Stretcher Lift.
- Ensure that the whirlpool inlet screen (23) is in place
- Turn the whirlpool outlet jet (21), anti-clockwise to full open position
- Fill bathtub automatically by pressing the on/off control (2) once.
- Set and monitor fill water temperature with mixer (3) and thermometer (8)

**Note:** Always physically check water temperature before and during bathing.

Read operating manual for patient handling.

#### Showering



- Lift the shower handle (11) on the right side on the panel with the white hose (10) and hold it over the tub, away from patient.
- Start the patient shower by pushing the on/off control (4) once on front of the control panel.
- Set the mixer (3) to desired temperature.
- Physically test water temperature and adjust the mixer to suit.
- Shower or rinse patient with shower spray handle.
- Stop showering by pushing the on/off control (4). Once and return handle to its bracket (12)

**Note:** Shower handle shall not be held in the bath water while bathing.

#### Whirlpool

- The water level shall cover the outlet jet (21).
- Adjust the jet (21) for desired flow direction: centred for perineal area; to either side for circular flow. Avoid placing the jet in a full upward position, as this could cause the water to splash.
- Make sure the whirlpool inlet screen (23) is not obstructed.
- With the patient in the bathtub, turn on the whirlpool unit by pushing the "ON-OFF" control button (5).
- To adjust water flow direction, turn the whirlpool jet (21).
- Increase or decrease whirlpool stream by turning the air control (22).
- Turn the whirlpool unit off by pushing the "ON-OFF" control button (5).

**Note:** This document is a complement to the Operating Manual, not a substitute.

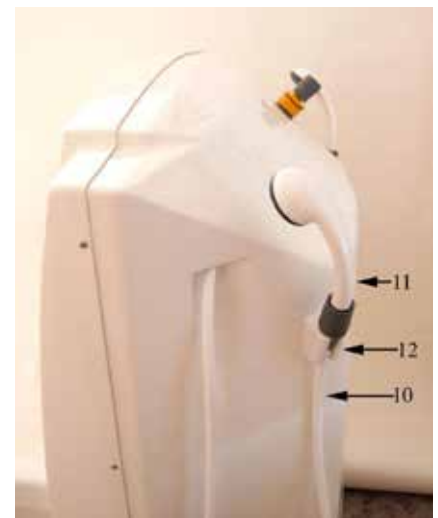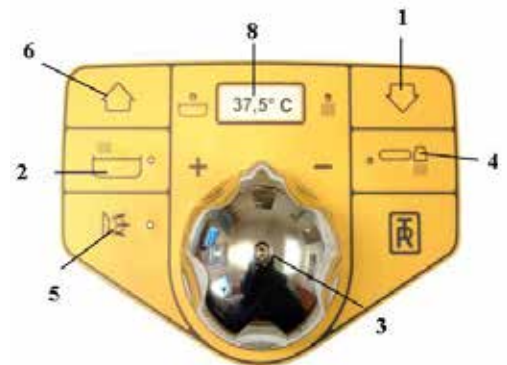

#### Head office


TR Equipment AB, Sweden  
Phone: +46 140 38 50 60  
info@trequipment.se  
www.trequipment.com

#### Group Subsidiaries

TR Group Inc, US, Phone: +1 (800) 752 6900  
TR Equipment Ltd, UK, Phone: +44 (0) 844 335 8386  
TRanspoRehaMed GmbH, Germany, Phone: +49 (0) 212 645850-0



TR EQUIPMENT

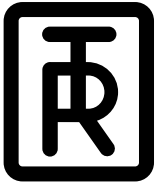

# TR 900A

## Height adjustable Autofill bathtubs


Hygiene Equipment  
Designed for People



### After bathing

#### Clean the equipment.


- Raise the tub to a comfortable height with the up control button (6). Remove the cleaning spray handle (15) with the yellow hose.
- Hold the handle over the bathtub while opening the shut-off safety valve (13) with its separate key. The cleaning detergent is automatically mixed with cold water to minimise inhalation.
- Briefly spray the cleaning mixture into the tub while checking the colour, smell and the foam of the mixture. Spray over all interior surfaces of the tub. Remove residue and debris with cloth or brush.
- To clean TR Chair Lift and Stretcher Lift, position them over the tub and spray them with cleaning mixture. The Chair and Stretcher have removable pads and the Stretcher has a head pillow which may be placed in the tub for cleaning.
- After cleaning, rinse with fresh water using the shower handle (11) with the white hose and hold it over the tub.
- Start the patient shower by pushing the on/off control (4) on panel.
- Set the mixer (3) to desired temperature.

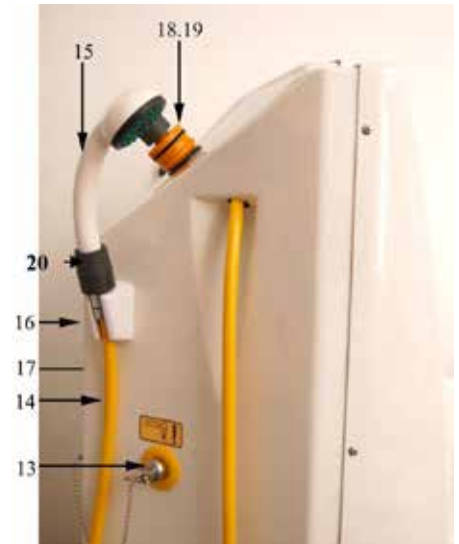

#### Clean the whirlpool system

A yellow inlet cleaning adapter (19) is stored on top of the detergent filling lid (18).

- Remove the whirlpool inlet safety screen (23) by turning the screen anti clockwise.
- Fit in the yellow adapter (19) to the whirlpool inlet (24).
- Remove cleaning shower handle (15) from yellow cleaning hose (14) by pressing the buttons on the quick coupling (20)
- Snap the hose onto the adapter and flush cleaning through the system by turning the cleaning "key" valve (13).
- Cleaning mixture will flush out through the whirlpool outlet (21) in the tub.
- When this is done, turn cleaning system off.
- Repeat procedure with patient shower handle (11) with the white hose to flush system with fresh water.
- Replace the yellow adapter on top of the detergent filling lid (18).

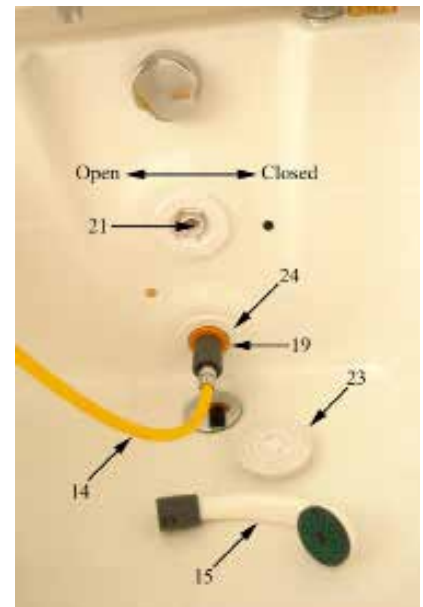

#### Detergent filling

If there is no cleaning liquid left, open the lid (18) on top of the panel and fill with maximum 1,5 litre (0,40 gallon) of concentrated detergent.

#### Study the Operating Manual for calibration of cleaning mixture and handling information.

#### Head office

TR Equipment AB, Sweden  
Phone: +46 140 38 50 60  
info@trequipment.se  
www.trequipment.com

#### Group Subsidiaries

TR Group Inc, US, Phone: +1 (800) 752 6900  
TR Equipment Ltd, UK, Phone: +44 (0) 844 335 8386  
TRanspoRehaMed GmbH, Germany, Phone: +49 (0) 212 645850-0



TR EQUIPMENT