



TR Classic

Tilt Table

*Hygiene Equipment
Designed for People*

TR Classic, Provides a Smooth and Effortless Lift

The TR Classic Tilt Table is designed to make testing, rehabilitation, positioning and transferring patients easier and safer. It provides a smooth and effortless lift from a horizontal to upright position.



CE



Registered at Medical Products
Agency, Sweden,
as a MD Class 1 product,
MDD 93/42EEC.


TR reserves the right to make
design changes and/or
improvements without notice.

The TR Classic Tilt Table is a logical choice for patients who have paralysis, are bedridden or wheelchair bound and could benefit from passive standing to improve circulation or reduction of lower extremity edema.

The TR Classic Tilt Table can be used in all environments to assist the caregiver with transfer and treatment needs. Standard equipment includes three safety belts for comfort and security, footboard, removable wire basket and a hand operated control.



TR EQUIPMENT

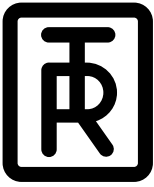

TR Classic Tilt Table

Hygiene Equipment
Designed for People

Features

- The TR Classic Tilt Table is electrically operated for smooth and silent operation. It effortlessly tilts from horizontal to a 90 degree upright position.
- A hand control is provided for convenient tilt adjustment and tilt indicator for precise monitoring of positioning.
- The wide footplate, adjustable and lockable in three angle positions, meets the varied needs of the patient population.
- The table comes standard with a back-up battery for lowering in case of a power failure. Wire storage basket is removable.
- The TR Classic Tilt Table easily maneuvers on four 125 mm (5") castors. Wheel fenders protect the lying surface when transporting through narrow doorways.



TR Classic has a wide foot plate, adjustable and lockable in three angle positions.



TR Classic Tilt Table is operated by an electrical actuator with battery back up for lowering in case of power failure.



TR Classic effortlessly tilts from horizontal to a 90 degree upright position and has a tilt indicator for monitoring of angle.

Specifications



| | |
|-------------------------------------|--|
| Cushion length - Standard | 1750 mm (69") |
| Cushion width - Standard | 650 mm (24 1/2") |
| Height (Horizontal position) | 600 mm (24") |
| Working Load Limit | 150 kg (330 lbs) |
| Power supply | 230 V, 50 Hz |
| Power failure lowering | Battery backup (9 V) |
| Frame | Powder coated steel |
| Castors | 125 mm (5") Two castors with brakes and two with straight steering |
| Control | Hand control |
| Cushion | Foam filled, vinyl cover |
| Models | Standard |
| Accessories | Activity Tray |

Head office
TR Equipment AB, Sweden
Phone: +46 140 38 50 60
info@trequipment.se
www.trequipment.com

TR Group Inc
130 North Madison Street
Rockford, IL61107, USA
Phone: +1 (800) 752 6900
Email: information@trequipment.com
www.trequipment.com



TR EQUIPMENT