



TR Equipment

Hygiene Equipment
Designed for People

Multi purpose trolleys brochure

CE



Registered at Medical Products
Agency, Sweden,
as a MD Class 1 product,
MDD 93/42EEC.

TR reserves the right to make
design changes and/or
improvements without notice.



Endless Possibilities

Showering
Transport
Dressing
Treatment



TR EQUIPMENT



Shower trolleys



TR Equipment offers the markets widest Shower trolley program with several options regarding sizes and Working Load Limits (WLL), both battery and hydraulic operated.

TR 2000

Hydraulic operated

Standard
Pediatric
Long

Working Load Limit: 150 kg (330 lbs)



The TR Shower Trolley is unique with its side mounted column, allowing the trolley stretcher to overlap the bed further for a safe and easy person transfer. The trolley assists and functions already at the bed side allowing the attendant to hands free access to the trolley controls thanks to foot operated height adjustment, castor brakes and straight steering.

The protective side rails are lowered when the stretcher is located over the bed and without any heavy lifting for the attendant, the person is easily, gently and safely positioned on the trolley mattress. The trolley side rails are then secured before transferring away from bed.

The side rail on the attendant side is hinged so that it will fold all the way down. This makes it easier for the attendant to handle the person during transfer. The trolley stretcher can be tilted in both directions for "Trendelenburg" position and has drains at both ends allowing the person to be positioned in either directions. There is no head or foot end specified.

The TR 2000/3000 can be delivered in three different lengths and the mattress and pillow are available in three different colours: Grey (Standard), Blue and Yellow (Option).




TR 3000

Battery operated

Standard
Pediatric
Long

Working Load Limit: 150 kg (330 lbs)



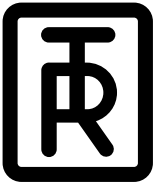

Bariatric shower trolleys



TR 4000 Atlas

Battery operated

Bariatric

Working Load Limit: 450 kg (1000 lbs)



The TR 4000/TR 4200 Multi Tasking Shower Trolley is designed for lots of power and loads of space for your most difficult patients. The trolley is used for shower, transport, as a changing or treatment trolley, the possibilities are endless.

The person can be transferred from bed to shower and back again in a safe and secure manner. Rechargeable battery powered dual actuators raise and lower the stretcher easily and quietly, with an easy to use hand control. The hand control operates the power tilt in both directions for improved drainage and tilts to a "Trendelenburg" position when desired.

An audible tone sounds when the stretcher is level. The side rails keep the person safely in place. Side rails swing down for transfer and lock into raised position for person security and a secure transport. The trolley has a foot controlled central locking castor brakes and straight steering. For a tall person the head end rail can, on TR 4000, be folded to a horizontal position for extra length. The stretcher has two positions, flat for transfer, tilt for transport shower and improved drainage.




TR 4200 Atlas Junior

Battery operated

Bariatric

Working Load Limit: 270 kg (530 lbs)



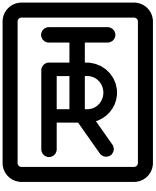

TR 2810 Shower panel



TR 2810 Shower panel has a thermoscopic mixing valve for the shower and on/off valves for the shower and the cleaning handles. The cleaning system runs on cold water to minimize inhalation when cleaning equipment. Partly built-in visible bottle with lid for cleaning detergent.

The TR Shower Trolleys make the showering and transferring easier for attendants. When integrated with the TR 2810 Shower Panel, TR presents a complete shower system.



Specifications

Models	Type	Length (mm/inch)	Width (mm/inch)	Height (mm/inch)	Working Load Limit (kg/lbs)
2000-S, 3000-S	Pediatric	1600 / 63	705 / 27 ^{3/4}	625-925 / 24 ^{1/2} - 36 ^{1/2}	150 / 330
2000-SH, 3000-SH	High chassis	1600 / 63	705 / 27 ^{3/4}	685-985 / 27-39	150 / 330
2000, 3000	Standard	1900 / 75	705 / 27 ^{3/4}	625-925 / 24 ^{1/2} - 36 ^{1/2}	150 / 330
2000-H, 3000-H	High chassis	1900 / 75	705 / 27 ^{3/4}	685-985 / 27-39	150 / 330
2000-L, 3000-L	Long	2100 / 83	705 / 27 ^{3/4}	625-925 / 24 ^{1/2} - 36 ^{1/2}	150 / 330
2000-LH, 3000-LH	High chassis	2100 / 83	705 / 27 ^{3/4}	685-985 / 27-39	150 / 330
4000 Atlas	Bariatric	2210 / 87	915 / 36	535-835 / 21-33	400 / 1000
4200 Atlas Junior	Bariatric	2150 / 85	740 / 29	535-835 / 21-33	240 / 530

All products are manufactured in Sweden

Distributor:

