



# TR 4200 Shower Trolleys

Hygiene Equipment  
Designed for People

## TR 4200, makes the showering and transferring easier for attendants

The TR 4200 Multi Tasking Bariatric Shower Trolley is designed for lots of power and loads of space for your most difficult resident needs. The TR 4200 has many application possibilities including showering, transporting and providing an efficient changing or treatment surface area. Potential impact wear areas are manufactured from powder coated stainless steel.



Registered at Medical Products Agency, Sweden, as a MD Class 1 product, MDD 93/42EEC.

TR reserves the right to make design changes and/or improvements without notice.



The resident can be transferred from bed to shower and back again in a safe and secure manner. The height of the base allows mobile sling lift transfers.

Rechargeable battery powered dual actuators raise and lower the stretcher easily and quietly with an easy to use hand control. The hand control operates the power tilt in both directions for improved drainage and tilts to a "Trendelenburg" position when desired.

An audible tone sounds when the stretcher is level.



**TR EQUIPMENT**



# TR 4200 Shower Trolleys

Hygiene Equipment  
Designed for People

## Features

- The stretcher is easily and quietly raised by two powerful actuators. The hand control operates the power tilt in both directions and raises and lowers the trolley.
- The stretcher can easily be tilted to the "Trendelenburg" position, to alleviate a person in shock. The height adjustable stretcher enables a comfortable working height position for the attendant. An audible tone sounds when the stretcher is level. When transporting through narrow doorways the side rails can be folded reducing the trolley stretcher width.
- The stretcher has rubber bumper wheels on all four corners to protect walls and doorways. Stretcher comes standard with a soft flame retardant mattress with drainage at the foot end and includes a drain hose and pillow.
- The trolley is made of durable powder coated steel, zinc plated steel that make cleaning simple and efficient.
- The mattress is easily removable.



Convenient single foot control for straight steering and brakes.

Side rails swing down or lock in horizontal position for transfer and upright position for transport.

## Specifications




<b>Length</b>	2150 mm (85")
<b>Width</b>	740 mm (29")
<b>Height</b>	535 - 835 mm (21" - 33")
<b>Working Load Limit</b>	270 kg / 600 lbs (Patient + Water)
<b>Brake bar and upper side rail/hinge</b>	Powder coated stainless steel
<b>Castors</b>	Brakes on all four castors and straight steering
<b>Control</b>	Hand control for height adjustment and tilting
<b>Mattress, Pillow</b>	Flame retardant material
<b>Battery</b>	Rechargeable, 24 V
<b>Accessories</b>	Backrest wedge pillow Backrest wedge pillow 15° Backrest wedge pillow 30°

## Accessories/Options

- All models can be delivered with grey or blue mattress.
- Fasteners/straps for attachment on mattresses to siderails.

## Dimensions (mm)

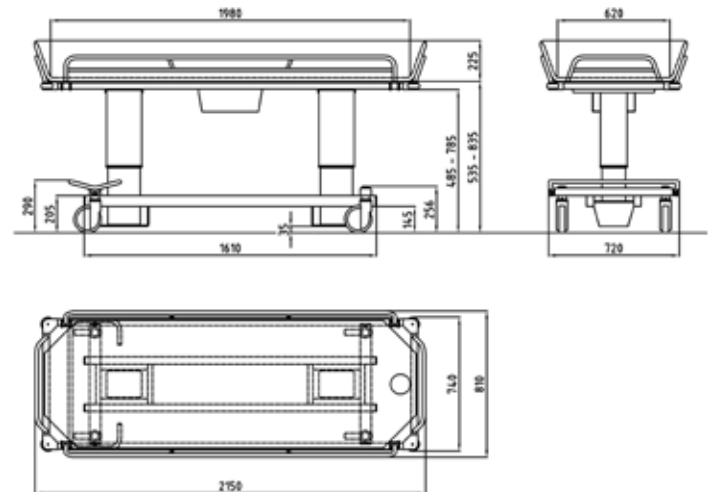

## Head office

TR Equipment AB, Sweden  
Phone: +46 140 38 50 60  
info@trequipment.se  
[www.trequipment.com](http://www.trequipment.com)

## Group Subsidiaries

TR Group Inc, US, Phone: +1 (800) 752 6900  
TR Equipment Ltd, UK, Phone: +44 (0) 844 335 8386  
TRanspoRehaMed GmbH, Germany, Phone: +49 (0) 212 645850-0



**TR EQUIPMENT**