


StaminaLift® TS6000 Transfer System

Power assisted manual handling solutions for your environment

The StaminaLift® TS6000 Transfer System offers proven technology designed to prevent work related injury by reducing the effort required to move hospital beds.

Features

- ✓ Universal compatibility
 - Designed for global Bariatric Bed compatibility
- ✓ Adjustable Jaws to suit Bariatric Bed
- ✓ Wheel base width: 550mm to 1010mm
 - Castor sizes: Ø125mm to Ø200mm (single and double)
 - Jaw bypass - releases bed in an emergency
- ✓ Brake bypass - allows freewheeling in an emergency
- ✓ 900kg Safe Working Load (SWL)
- ✓ Adjustable handle to suit operator height
- ✓ Works on all indoor surfaces, non-marking wheels
- ✓ Battery charger suits 120V and 240V AC supply
- ✓ Unique patented technology
- ✓ Dual drive 360-degree movement
- ✓ Keyless control with safety lockout
- ✓ One-person operation
 - Hands free connect and lock
 - Fingertip control
 - No muscular strain
- ✓ Quiet Operation: 65dB
- ✓ Automatic braking

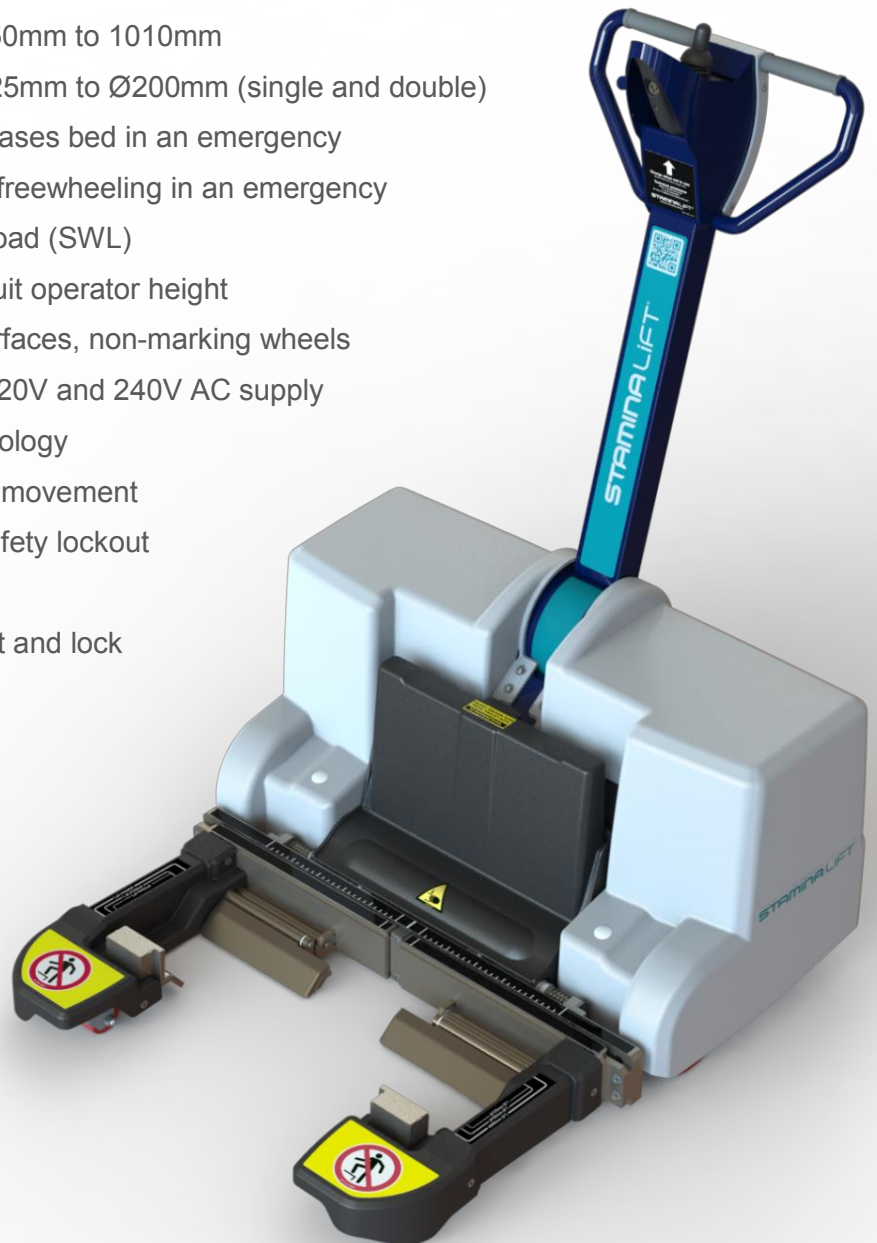



STAMINALIFT

StaminaLift® TS6000 Transfer System

Technical Specifications

Max load capacity	450 kg lift and 900 kg push or pull 993 lbs lift and 1985 lbs push or pull
Weight	200 kg/441 lbs
Dimensions (L x W x H)	In Use: 820mm x 900mm x 1225mm - (33in. x 36in. x 48in.) Packaged: 850mm x 940mm x 480mm - (34 in. x 37 in. x 19 in.)
Batteries	Sealed maintenance free GEL batteries.
Recharge	Nightly recharge recommended using supplied automatic charger
Warning device	Horn (warning device) on control unit
Security	Electronic locking function available, refer to operators manual
Wheels	Non-marking, suitable for indoor use on a variety of surfaces
Controls	Programmable, ergonomic joystick for control of direction, speed and jaw movement
Safety features	Automatic brakes / Emergency shut-down button Emergency Jaw release – allows override of lowering mechanism Emergency Brake release – allows override of brake mechanism Safety features to prevent unauthorized use
Motor	Two variable drive DC electric motors
Servicing	Semi-Annual



TR Group Inc
130 North Madison Street
Rockford, IL 61107, USA
Phone: +1 (800) 752 6900
Email: information@trequipment.com
www.trequipment.com

STAMINALIFT®

www.staminalift.com

StaminaLift International Limited
Tel: +61 8 8262 7137
Fax: +61 8 8359 5745
Email: info@staminalift.com