



# TR Classic

## Tilt table

Hygiene Equipment  
Designed for People

### TR Classic, Provides a Smooth and Effortless Lift

The TR Classic Tilt Table is designed to make testing, rehabilitation, positioning and transferring patients easier and safer. It provides a smooth and effortless lift from a horizontal to upright position.

CE



Registered at Medical Products  
Agency, Sweden,  
as a MD Class 1 product,  
MDD 93/42EEC.

TR reserves the right to make de-  
sign changes and/or improvements  
without notice.




The TR Classic Tilt Table is a logical choice for patients who have paralysis, are bedridden or wheelchair bound and could benefit from passive standing to improve circulation or reduction of lower extremity edema.

The TR Classic Tilt Table can be used in all environments to assist the caregiver with transfer and treatment needs. Standard equipment includes three safety belts for comfort and security, footboard, removable wire basket and a hand operated control.



**TR EQUIPMENT**

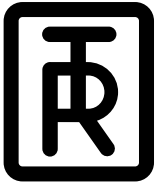

# TR Classic

## Tilt table

Hygiene Equipment  
Designed for People

### Features

- The TR Classic Tilt Table is electrically operated for smooth and silent operation. It effortlessly tilts from horizontal to a 90 degree upright position.
- A hand control is provided for convenient tilt adjustment and tilt indicator for precise monitoring of positioning.
- The wide footplate, adjustable and lockable in three angle positions, meets the varied needs of the patient population.
- The table comes standard with a back-up battery for lowering in case of a power failure. Wire storage basket is removable.
- The TR Classic Tilt Table easily maneuvers on four 125 mm (5") castors. Wheel fenders protect the lying surface when transporting through narrow doorways.



TR Classic has a wide foot plate, adjustable and lockable in three angle positions.



TR Classic Tilt Table is operated by an electrical actuator with battery back up for lowering in case of power failure.



TR Classic effortlessly tilts from horizontal to a 90 degree upright position and has a tilt indicator for monitoring of angle.

TR Classic effortlessly tilts from horizontal to a 90 degree upright position and has a tilt indicator for monitoring of angle.

### Specifications



<b>Cushion length - Standard</b>	1750 mm (69")
<b>Cushion width - Standard</b>	650 mm (24 1/2")
<b>Height (Horizontal position)</b>	600 mm (24")
<b>Working Load Limit</b>	150 kg (330 lbs)
<b>Power supply</b>	230 V, 50 Hz
<b>Power failure lowering</b>	Battery backup (9 V)
<b>Frame</b>	Powder coated steel
<b>Castors</b>	125 mm (5") Two castors with brakes and two with straight steering
<b>Control</b>	Hand control
<b>Cushion</b>	Foam filled, vinyl cover
<b>Models</b>	Standard
<b>Accessories</b>	Activity Tray

### Head office

TR Equipment AB, Sweden  
Phone: +46 140 38 50 60  
info@trequipment.se  
[www.trequipment.com](http://www.trequipment.com)

### Group Subsidiaries

TR Group Inc, US, Phone: +1 (800) 752 6900  
TR Equipment Ltd, UK, Phone: +44 (0) 844 335 8386  
TRanspoRehaMed GmbH, Germany, Phone: +49 (0) 212 645850-0



TR EQUIPMENT